

AUTOMATIC TUBE AND VIAL LABELING WITHOUT THE FRONT-END BOTTLENECK

FLEX TVF (FLEXIBLE TUBE AND VIAL FEEDER)

The new Flex TVF makes your PR100i + TubePro Accessory even more productive—it eliminates the process of manually presenting tubes and vials one-by-one for labeling.


Printing and applying labels automatically to tubes and vials is a huge step forward in laboratory automation. But the tube still needs to be positioned correctly in order to be labeled, which means manual labor. Now that step can be automated as well, with the new Flex TVF.

The Flex TVF is a highly automated bulk feeder unlike any other. It will automatically separate individual tubes from bulk storage and ensure that each tube exits the system in exactly the right orientation for labeling. Orienting and dispensing more quietly than any other feeder, tubes and vials of various shapes and sizes can be handled- including those with rubber rings or gaskets around the neck.

Extend your lab automation beyond what you ever thought possible.

Features:

- › Handles tubes and vials of various sizes and shapes
- › Efficient processing at more than twice the speed of manual labeling
- › Replaces antiquated vibratory bowls because:
 - TVF is much quieter
 - Unlike vibratory bowls that lose some product during operation, the TVF handles 100% of the product loaded
 - Changeover from one vial size to another is easy
 - Vials with a rubber ring or gasket can be handled, which would clog a standard vibratory bowl
 - TVF can handle straight tubes with no lips, with or without caps; vibratory bowls need some sort of “lip” on them to feed through the channel in the bowl
 - TVF is smart enough to stop processing if, for example, an uncapped vial is inadvertently included in a batch of capped vials

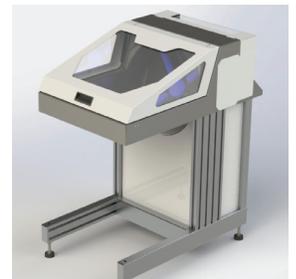

SAVE LIVES, SAVE TIME, SAVE MONEY.

There's a lot of pressure put on the labels you use. If they don't perform well when it comes to readability, trackability and just hanging on, there's more at stake than just a missing label. With adhesives and materials that can withstand harsh environments, and barcode symbols that have been checked and double-checked, Computype's got you covered with products that perform.

LOOK TO US FOR ALL YOUR CRITICAL LABELING SOLUTIONS.

We work closely with you to provide the right labels to meet all your requirements. And not just the right labels, but everything for your labeling needs. From labels to printers, and scanners to label replication systems, Computype offers this and everything in between.

WE'D LOVE TO HEAR FROM YOU.

Computype
2285 West County Road C
St. Paul, MN 55113-2567
Phone: 800-328-0852 | 651-633-0633
Fax: 651-633-5580
Email: sales@computype.com
Twitter: [@computype](https://twitter.com/computype)

www.computype.com

Computype
Kingswood Business Park, Connaught Road
Kingston Upon Hull
HU7 3AP
United Kingdom
Phone: [44] 1482-835366
Fax: [44] 1482-831152 421-606
Email: contactus@computype.com

LS0923.r2.0417



Computype™